



Veinlite®



PEDI2

For Pediatrics



14 bright light emitting diodes (LEDs)

4 deep red LEDs

10 high-contrast orange LEDs



Ideal for the pediatric patients, the Veinlite PEDI2® is specifically designed for infants through adolescents. Powered by two AAA batteries, this smallest and most compact version of the Veinlite LED product line has two switches operating 14 dual colored LEDs for maximum contrast. The PEDI2 has 3 different levels of brightness for variable light output.

FEATURES

- Small compacted size for pediatric patients
- Color selection
- 3 levels of brightness
- Auto shut-off timer

KIT INCLUDES:


- Veinlite PEDI2®
- 25 disposable covers
- 2 AAA batteries

NEO

For Neonates



Orange Light
for vein viewing

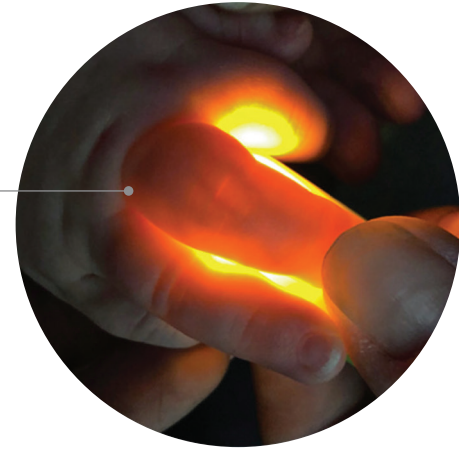

Veinlite® NEO is a revolutionary transillumination device designed to detect both veins and arteries in neonates. Veinlite NEO's new form factor makes it easy, convenient and readily available for assisting with venous and arterial access in neonatal pediatric patients. This device features 3 different colored lights: green for arteries, orange for vein imaging, and a white exam light.

FEATURES

- Only device on the market for arterial sticks
- New, compact form factor
- 3 light colors and 3 levels of brightness
- Through-the-body transillumination
- Auto shut-off timer

KIT INCLUDES:


- Veinlite NEO®
- Lanyard-ready plastic cap
- AAA 1.5v battery
- 50 disposable plastic covers

R

For Research



Orange Light
for viewing
arteries; shown
with rat tail adapter

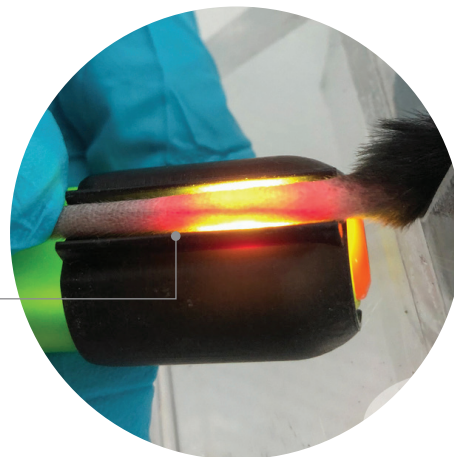

Veinlite R® is a revolutionary transillumination device designed to detect both veins and arteries in mice and rats. Veinlite R's new form factor makes it easy, convenient and readily available for assisting with venous and arterial access. The included tail adapters make stabilizing and securing the tail during access effortless. This device features 3 different colored lights: green for arteries, orange for vein imaging, and a white exam light.

FEATURES

- Only device on the market for arterial sticks
- New, compact form factor
- 3 light colors and 3 levels of brightness
- Through-the-body transillumination
- Auto shut-off timer

KIT INCLUDES:

- Veinlite R®
- Lanyard-ready plastic cap
- AAA 1.5v battery
- 50 disposable plastic covers

Included Separate bag:

- Large Rat Tail Adapter
- Small Mouse Tail Adapter



Product Comparison



Product Name	VEMS-PRO	VLED+	VLEDX	VPEDI2	VNEO	VLR
Number of LEDs	24	28	32	14	3	3
Color of LEDs	Orange (12) + Red (4) + White (8)	Orange (22) + Red (6)	Orange (24) + Red (8)	Orange (10) + Red (4)	Orange (1) + Green (1) + White (1)	Orange (1) + Green (1) + White (1)
Color Selection	No	Yes	Yes	Yes	Yes	Yes
Ring Diameter (mm)	21	27	31	18	N/A	N/A
Battery	2xAA Alkaline	Rechargeable Lithium Ion	Rechargeable Lithium Ion	2 AAA 1.5v	AAA 1.5v	AAA 1.5v
Brightness Control	No	Yes	No	Yes	Yes	Yes
Dimensions	115 x 55 x 21	99 x 57 x 21	102 x 65 x 21	85 x 40 x 17	96.1x19x16.5	96.1x19x16.5
Weight, Inc. Battery (g)	108	77	83	64	29	29

**Testing performed on low power setting*

To order your Veinlite visit veinlite.com or call (281) 240-3111.

Patent Number US 7,874,698 B2 | © Registered Trademark of TransLite, LLC Designed and manufactured in the USA.